



### Healthcare Solutions



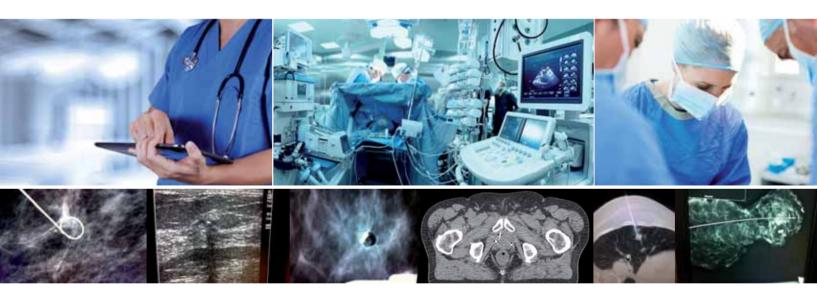

**EXPORT CATALOGUE** 

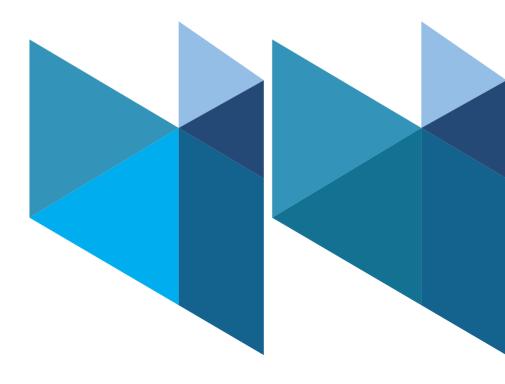





### **COMPACT** SOFT TISSUE BIOPSY

#### Programmable Automatic Biopsy System



- One-handed design simplifies loading and operation whilst allowing for convenient manipulation of the ultrasound transducer
- Automatic mode: releases cannula and stylet in rapid sequence and captures tissue sample with the push of a button
- Semiautomatic: releases the stylet allowing time for tissue to settle in notch (ideal for denser tissue or when there is a desire to verify notch position), and activates the cutting cannula
- 20mm sample notch and thin-wall cannula for larger core
- Colour coded for guick gauge size identification
- Depth markings etched on cannula for precise positioning
- The Coaxial outer needle permits multiple biopsies (US and CT guided)
- · Commonly used for the following biopsy procedures: Breast, Kidney, and Lung



One-handed design simplifies loading and operation whilst allowing for convenient manipulation of the ultrasound transducer

COMPACT Qty/Box 10pcs					
Add C for Coaxial	Add <b>C</b> for Coaxial set ( VCTC )				
Add 1 for Trocar t	ip ( VCT1 )				
Add <b>A</b> for Prostate	e kit ( VCTA )				
CODE	SIZE	COAXIAL KIT			
• VCT1210	12Gx10cm	VCTC1210			
● VCT1410	14Gx10cm	VCTC1410			
● VCT1415	14Gx15cm	VCTC1415			
VCT1615	16Gx15cm	VCTC1615			
<ul><li>VCT1810</li></ul>	18Gx10cm	VCTC1810			
<ul><li>VCT1815</li></ul>	18Gx15cm	VCTC1815			
<ul><li>VCT1820</li></ul>	18Gx20cm	VCTC1820			
<ul><li>VCT1825</li></ul>	18Gx25cm	VCTC1825			
<ul><li>VCT2010</li></ul>	20Gx10cm	VCTC2010			
<ul><li>VCT2015</li></ul>	20Gx15cm	VCTC2015			



#### Automatic Biopsy Device



- One-handed design simplifies loading and operation whilst allowing for convenient manipulation of the ultrasound transducer
- Equipped with a visual and mechanical safety-system
- 20mm sample notch and thin-wall cannula for larger core
- Two triggers permit ergonomic handling for right and left-handed application, regardless of the position of the extraction location
- Colour coded for quick gauge size identification
- Depth markings etched on cannula for precise positioning
- The Coaxial outer needle permits multiple biopsies (US and CT guided)
- Commonly used for the following biopsy procedures: Breast, Kidney, and Prostate



Mechanical safety-system



One-handed design

COREX Qty/Box 10 <i>pcs</i>					
Add C for Coaxi	Add <b>C</b> for Coaxial set ( VBSC )				
Add 1 for Trocar	tip ( VBS1 )				
Add <b>A</b> for Prosta	nte kit ( VBSA )				
CODE	SIZE	COAXIAL KIT			
● VBS1410	12Gx10cm	VBSC1210			
● VBS1415	14Gx15cm	VBSC1415			
VBS1610	16Gx10cm	VBSC1610			
○ VBS1615	16Gx15cm	VBSC1615			
● VBS1810	18Gx10cm	VBSC1810			
● VBS1815	18Gx15cm	VBSC1815			
● VBS1820	18Gx20cm	VBSC1820			
● VBS1825	18Gx25cm	VBSC1825			
<ul><li>VBS2010</li></ul>	20Gx10cm	VBSC2010			
<ul><li>VBS2015</li></ul>	20Gx15cm	VBSC2015			



### CT-CORE SOFT TISSUE BIOPSY

#### Single Action Biopsy Device



- Lightweight and compact to facilitate one-handed operation during procedures
- Precise, rapid-firing mechanism for an intact core tissue sample
- Specimen notch fully exposed when stylet is advanced for ease of sample removal
- The puncture depth can be adjusted between 15 to 22 mm
- Beveled stylet for easy penetration into the lesion with less trauma to surrounding tissues and organs
- Colour coded for quick gauge size identification
- The 14 gauge version features a Trocar point for ease of insertion into dense tissue
- Echogenic tip for enhanced visualization under ultrasound guidance
- Depth markings etched on cannula for precise positioning
- The Coaxial outer needle permits multiple biopsies (US and CT guided)
- Commonly used for the following biopsy procedures: Kidney, Liver, and Lung



One-handed design



Detachable cannula (VVD0A)

CT-CORE Qty/Box 10 <i>pcs</i>					
Add C for Coaxia	Add <b>C</b> for Coaxial set ( VVD0C )				
Add 1 for Trocar	tip (VVD01)				
Add <b>A</b> for Detac	hable cannula ( VVI	D0A )			
CODE	SIZE	COAXIAL KIT			
• VVD01410	14Gx10cm	VVD0C1410			
● VVD01415	14Gx15cm	VVD0C1415			
VVD01610	16Gx10cm	VVD0C1610			
VVD01615	16Gx15cm	VVD0C1615			
<ul><li>VVD01810</li></ul>	18Gx10cm	VVD0C1810			
<ul><li>VVD01815</li></ul>	18Gx15cm	VVD0C1815			
<ul><li>VVD01820</li></ul>	18Gx20cm	VVD0C1820			
<ul><li>VVD01825</li></ul>	18Gx25cm	VVD0C1825			
<ul><li>VVD02010</li></ul>	20Gx10cm	VVD0C2010			
<ul><li>VVD02015</li></ul>	20Gx15cm	VVD0C2015			



#### Reusable Biopsy Instrument



- Equipped with a visual and mechanical safety-system
- Two triggers permit ergonomic handling for right and left-handed application, regardless of the position of the extraction location
- The puncture depth can be adjusted between 16 to 22 mm
- Compatible exclusively with Vigeo Autocut Biopsy Needle
- · Loaded status indicator





Loaded status indicator



**AUTOCUT V** 

Specially designed to use with VTEK®-VVTEK®

V-TEK				
Add $oldsymbol{V}$ for Single trigger version ( VV-TEK )				
CODE	THROW	SAFETY		
V-TEK	16-22 mm	YES		
VV-TEK	22 mm	NO		
AUTOCUT Qty/Box	( 10 <i>pcs</i>			
Add C for Coax	ial set ( VVGC )			
Add 1 for Troca	r tip ( VVG1 )			
Add <b>A</b> for Prosta	ate kit (VVGA)			
CODE	SIZE	COAXIAL KIT		
<ul><li>VVG1210</li></ul>	12Gx10cm	VVGC1210		
● VVG1410	14Gx10cm	VVGC1410		
<ul><li>VVG1413</li></ul>	14Gx13cm	VVGC1413		
VVG1610	16Gx10cm	VVGC1610		
VVG1615	16Gx15cm	VVGC1615		
<ul><li>VVG1810</li></ul>	18Gx10cm	VVGC1810		
<ul><li>VVG1815</li></ul>	18Gx15cm	VVGC1815		
<ul><li>VVG1820</li></ul>	18Gx20cm	VVGC1820		
<ul><li>VVG1825</li></ul>	18Gx25cm	VVGC1825		
<ul><li>VVG2010</li></ul>	20Gx10cm	VVGC2010		
<ul><li>VVG2015</li></ul>	20Gx15cm	VVGC2015		



### **AUTOCUT** SOFT TISSUE BIOPSY

#### Refill Needle for Reusable Biopsy Instrument



AUTOCUT V Specially designed to use with VTEK®- VVTEK® Vigeo Device - MAGNUM® Bard Device



AUTOCUT P Specially designed to use with VVVTEK® Vigeo Device - PRO-MAG® ULTRA Argon Device



AUTOCUT E Specially designed to use with PRO-MAG® 2.2 Argon Device - BIOPTY® Bard Device

- Single patient use core biopsy needle designed for use with Reusable Biopsy Instrument
- Available in several needle gauge sizes and lengths
- Centimeter depth markings assist in needle placement
- Attached removable spacer maintains the distance between the two needle hubs promoting ease ofhandling and simplifying targetinge
- The Coaxial outer needle permits multiple biopsies (US and CT guided)
- Colour coded for quick gauge size identification
- Intended for use in obtaining biopsies from soft tissues such as Breast, Kidney, and Prostate

	AUTOCUT o	/D 40			
/	AUTOCUT Qty/Box 10pcs				
	Add <b>C</b> for C	oaxial set			
	Add 1 for Tr	ocar tip			
	Add <b>A</b> for Pr	ostate kit			
	CODE			SIZE	
	<ul><li>VVG1210</li></ul>	VVGP1210	VVE1210	12Gx10cm	
	● VVG1410	VVGP1410	VVE1410	14Gx10cm	
	● VVG1415	VVGP1415	VVE1415	14Gx15cm	
	VVG1615	VVGP1615	VVE1615	16Gx15cm	
	<ul><li>VVG1810</li></ul>	VVGP1810	VVE1810	18Gx10cm	
	<ul><li>VVG1815</li></ul>	VVGP1815	VVE1815	18Gx15cm	
	<ul><li>VVG1820</li></ul>	VVGP1820	VVE1820	18Gx20cm	
	<ul><li>VVG1825</li></ul>	VVGP1825	VVE1825	18Gx25cm	
	<ul><li>VVG2010</li></ul>	VVGP2010	VVE2010	20Gx10cm	
	<ul><li>VVG2015</li></ul>	VVGP2015	VVE2015	20Gx15cm	



# **HEPASHOT** SOFT TISSUE BIOPSY

#### One-handed Menghini Aspiration Device



- Single-handed safe and fast operation
- Pyramidal Tip allows for easy insertion
- Menghini design cannula for cylindrical and large sampling
- Internal echogenic marker to prevent cancerous cells spreading
- Intended for use in Liver, Lymph Nodes, and Thyroid
- Colour coded for quick gauge size identification
- Centimeter depth markings assist in needle placement
- The Coaxial outer needle permits multiple biopsies (US and CT guided)



One-handed design simplifies loading and operation whilst allowing for convenient manipulation of the ultrasound transducer

HEPASHOT Qty/Box 10pcs				
Add <b>C</b> for Coaxial set ( VVC0C )				
SIZE	COAXIAL KIT			
16Gx15cm	VVC0C1615			
16Gx20cm	VVC0C1620			
18Gx10cm	VVC0C1810			
18Gx15cm	VVC0C1815			
18Gx20cm	VVC0C1820			
20Gx05cm	VVC0C2005			
20Gx10cm	VVC0C2010			
20Gx15cm	VVC0C2015			
21Gx05cm	VVC0C2105			
21Gx10cm	VVC0C2110			
	SIZE  16Gx15cm 16Gx20cm 18Gx10cm 18Gx15cm 18Gx20cm 20Gx05cm 20Gx10cm 20Gx15cm 21Gx05cm			



### **(€** 0476

# HEPATIC SOFT TISSUE BIOPSY

#### Menghini Aspiration Device



- Pyramidal Tip allows for easy insertion
- Menghini design cannula for cylindrical and large sampling
- Internal echogenic marker to prevent cancerous cells spreading
- Intended for use in Liver, Lymph Nodes, and Thyroid
- · Centimeter depth markings assist in needle placement

HEPATIC Qty/Box 10	Opcs
CODE	SIZE
<ul><li>VVC1815</li></ul>	18Gx15cm
<ul><li>VVC1820</li></ul>	18Gx20cm
<ul><li>VVC2015</li></ul>	20Gx15cm

On request Vigeo is available to provide needles with different diameter and length

**€** 0476

### HEPASET SOFT TISSUE BIOPSY

Menghini Aspiration Device



- Biopsy needle, 90 mm
- Special Menghini bevel for excellent tissue sampling
- Thin-walled needle for minimal puncture trauma
- Aspiration Syringe 10 ml Injection Needle Scalpel

HEPASET Qty/Box 10pcs		
CODE	SIZE	
VVA1609	16Gx9cm	
<ul><li>VVA1709</li></ul>	17Gx9cm	
<ul><li>VVA1809</li></ul>	18Gx9cm	



# **EASYCUT** SOFT TISSUE BIOPSY

### Manual Capture Biopsy Device



- Manual advancement of cannula and stylet
- Removable stylet allows multiple sampling capability
- Thin-wall needle yields larger core
- Centimeter depth markings assist in needle placement




Removable stylet


EASYCUT Qty/Box 10pcs			
Add C for Coax	ial set ( VVMC )		
CODE	SIZE	COAXIAL KIT	
<ul><li>VVM1410</li></ul>	14Gx10cm	VVMC1410	
● VVM1415	14Gx15cm	VVMC1415	
VVM1610	16Gx10cm	VVMC1610	
VVM1615	16Gx15cm	VVMC1615	
<ul><li>VVM1810</li></ul>	18Gx10cm	VVMC1810	
<ul><li>VVM1815</li></ul>	18Gx15cm	VVMC1815	
<ul><li>VVM1820</li></ul>	18Gx20cm	VVMC1820	
<ul><li>VVM1825</li></ul>	18Gx25cm	VVMC1825	
<ul><li>VVM2010</li></ul>	20Gx10cm	VVMC2010	
<ul><li>VVM2015</li></ul>	20Gx15cm	VVMC2015	



### Vacuum Assisted Ultrasound Guided Breast Biopsy System



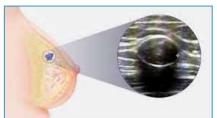

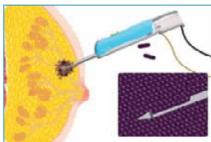

- Developed to overcome the weaknesses with complete biopsies
- Large amounts of samples can be obtained through a single insertion, thus allowing for a more accurate diagnosis
- Evaluated in clinical trials as a therapeutic tool, with an expanding role in the surgical treatment for excision of benign breast lesions such as fibroadenomas, nipple discharge, abscesses, and gynaecomastia
- An intuitive interface provides the operator easy control
- The aspirated tissue is automatically discharged by vacuum suction
- The inner blade has 4 step adjustments 5mm ,10mm ,15mm, and 20mm



X-Ray



Ultrasound

BEXCORE SYSTEM			
System body -	System body - Driver - Foot switch - Suction bottle		
CODE			
BXS100			
BEXCORE PRO	BE		
CODE	SIZE		
BXC145	8Gx109mm		
BXC140	10Gx109mm		
BEXCORE TISS	SUE MARKER		
CODE	SIZE	MARKER	
040239	363,7mm	2x4 mm	



### **BIOMARC** WOMEN'S HEALTH

#### Breast Tissue Marker Device







- Bar-bell shaped pyrolytic carbon coated ceramic with exceptionally biocompatibility with clear imaging characteristics
- Unique shapes designed to prevent migration
- No allergy response or hypersensitivity reported
- Designed for safe, reliable, and precise placement during core biopsies
- Permanent and clear visibility with mammography, ultrasound and MRI
- No interference with MRI visualization, lower level of artifact than metal markers
- No interference with Tomosynthesis visualization
- Available in three distinct size configurations
- Available for vacuum assisted and core biopsy needle devices

#### WHEN TO USE

- Marking after removal of tumor
- Marking microcalcification area after breast vacuum assisted biopsy (VAB)
- Marking of suspicious breast lesions after critical core biopsy
- Orientation for planning radiation therapy



Ultrasound



Mammography

BIOMAR	Qty/Box 10 <i>pcs</i>		
CODE	SIZE	MARKER	DESCRIPTION
040220	8Gx165,3mm	2x4mm	Flexible Applicator, 8G Mammotome®
040203	9Gx165,3mm	2x4mm	Flexible Applicator, 11G Mammotome®
040204	9 G x 1 2 0 m m	2x4mm	Flexible Applicator, 9G Eviva® / 10G Spirotome®
040228	9Gx363,7mm	2x4mm	Flexible Applicator, 9G Atec®
040239	9Gx363,7mm	2x4mm	Flexible Applicator, 10G SenoRX® / 10G Bexcore®
040306	13Gx117,7mm	2x4mm	Rigid Applicator, 13G
040320	14Gx117,7mm	1x4mm	Rigid Applicator, 14G
040300	14Gx117,7mm	1x3mm	Rigid Applicator, 14G
040310	17Gx117,7mm	0,9x5mm	Rigid Applicator, 17G




## RELOCK SOFT WOMEN'S HEALTH

Soft Repositionable Breast Localization Device

PATENT PENDING



- Preoperative marking of local lesions
- One-handed design for a convenient manipulation of the ultrasound transducer
- Retractable hook with unique fexibility of its Polyamid guiding thread
- Indicator on wing shows direction of hook
- Echogenic needle tip for ultrasound use
- Hook release at the needle tip for accurate targeting
- Centimeter depth markings assist in needle placement
- Patented stopper for partial hook release prior to confirming position
- 9mm Loop hook

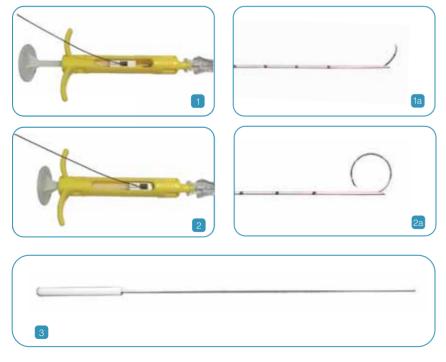

1. Half exposure - 2. Full exposure - 3. Stiffener cannula



Mammography

RELOCK SOFT Qty/Box 10pcs		
CODE	SIZE	
<ul><li>VVRPS1905</li></ul>	19Gx5cm	
<ul><li>VVRPS1907</li></ul>	19Gx7cm	
<ul><li>VVRPS1910</li></ul>	19Gx10cm	



# RELOCK WOMEN'S HEALTH

### Repositionable Breast Localization Device



- Preoperative marking of local lesions
- Retractable hook with unique fexibility
- Indicator shows direction of hook
- Echogenic needle tip for ultrasound use
- Hook release at the needle tip for accurate targeting
- Centimeter depth markings assist in needle placement
- Fixation piece for preoperative fixing of the wire
- 9mm Loop hook



Mammography

RELOCK Qty/Box 10pcs		
CODE	SIZE	
<ul><li>VVR1905</li></ul>	19Gx5cm	
<ul><li>VVR1907</li></ul>	19Gx7cm	
<ul><li>VVR1910</li></ul>	19Gx10cm	
<ul><li>VVR2007</li></ul>	20Gx07cm	
<ul><li>VVR2010</li></ul>	20Gx10cm	
<ul><li>VVR2015</li></ul>	20Gx15cm	




### Repositionable Breast Localization Device



- Preoperative marking of local lesions
- Retractable hook with unique fexibility
- Indicator shows direction of hook
- Echogenic needle tip for ultrasound use
- Hook release at the needle tip for accurate targeting
- Centimeter depth markings assist in needle placement
- Fixation piece for preoperative fixing of the wire
- Double 7mm Loop hook



Half exposure

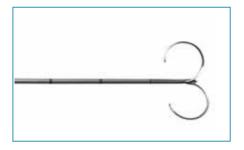

Full exposure

RELOCK DUO Qty/Box 10pcs		
CODE	SIZE	
<ul><li>VVNSP1905</li></ul>	19Gx5cm	
<ul><li>VVNSP1907</li></ul>	19Gx7cm	
<ul><li>VVNSP1910</li></ul>	19Gx10cm	



# RELOCK PREMIUM WOMEN'S HEALTH

Repositionable One-Handed Breast Localization Device

PATENT PENDING



- Preoperative marking of local lesions
- One-handed design for a convenient manipulation of the ultrasound transducer
- Removable handle for mammography procedure
- Retractable hook with unique fexibility
- Indicator shows direction of hook
- Echogenic needle tip for ultrasound use
- Hook release at the needle tip for accurate targeting
- · Centimeter depth markings assist in needle placement
- Patented stopper for partial hook release prior to confirming position
- 9mm Loop hook













RELOCK PREMIUM Qty/Box 10pcs

SIZE

19Gx5cm

19Gx7cm


19Gx10cm

CODE

VVRP1905

VVRP1907

VVRP1910

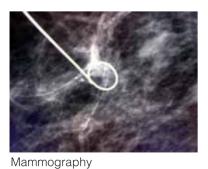




On request Vigeo is available to provide needles with different diameter and length Removable handle





- Preoperative marking of local lesions
- One-handed design for a convenient manipulation of the ultrasound transducer
- Removable handle for mammography procedure
- Retractable hook with unique fexibility
- Indicator shows direction of hook
- Echogenic needle tip for ultrasound use
- Hook release at the needle tip for accurate targeting
- Centimeter depth markings assist in needle placement
- Patented stopper for partial hook release prior to confirming position
- Double 7mm Loop hook



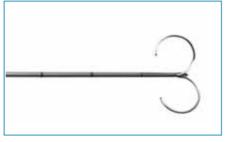

Half exposure

Full exposure




Removable handle

RELOCK PREMIUM DUO Qty/Box 10pcs		
CODE	SIZE	
VVNSPX1905	19Gx5cm	
<ul><li>VVNSPX1907</li></ul>	19Gx7cm	
<ul><li>VVNSPX1910</li></ul>	19Gx10cm	



## SURGIMARC WOMEN'S HEALTH

#### **Breast Localization Device**

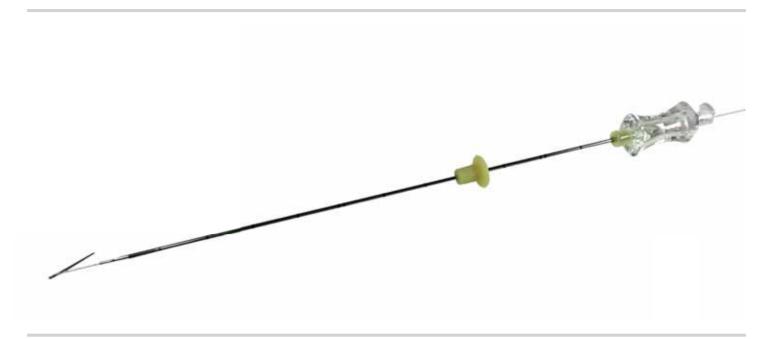

- Preoperative marking of local lesions
- Puncture cannula with ultrasonic tip with guidance wire in place
- Available Localization wire with reinforcement
- Transparent cannula connector
- The X-shaped hook provides excellent fixation and is intended to limit migration
- The 3 cm long stiff portion of the shaft is easily palpated during surgery
- Beads placed at 1cm intervals produce a rapid visual reference to the lesion
- Blunt Tip Stiffening Cannula available



SURGIMARC Qty/Box 20pcs				
Add <b>T</b> fo	r Stiffened Wi	re ( VVNT )		
Add <b>M</b> fo	or Wire Depth	Distal Markin	g (VVNM)	
Add <b>X</b> fo	r Double Anc	horing Hook (	VVNX)	
CODE				SIZE
<ul><li>VVN1807</li></ul>	VVNT1807	VVNM1807	VVNX1807	18Gx7cm
VVN1810	VVNT1810	VVNM1810	VVNX1810	18Gx10cm
<ul><li>VVN2007</li></ul>	VVNT2007	VVNM2007	VVNX2007	20Gx7cm
<ul><li>VVN2010</li></ul>	VVNT2010	VVNM2010	VVNX2010	20Gx10cm
<ul><li>VVN2107</li></ul>	VVNT2107	VVNM2107	VVNX2107	21Gx7cm
<ul><li>VVN2110</li></ul>	VVNT2110	VVNM2110	VVNX2110	21Gx10cm



## SURGIMARC SLIM WOMEN'S HEALTH

#### One-Handed Breast Localization Device



- Preoperative marking of local lesions
- One-handed design for a convenient manipulation of the ultrasound transducer
- Puncture cannula with ultrasonic tip with guidance wire in place
- Available Localization wire with reinforcement
- The X-shaped hook provides excellent fixation and is intended to limit migration
- The 3 cm long stiff portion of the shaft is easily palpated during surgery
- Beads placed at 1cm intervals produce a rapid visual reference to the lesion

• Blunt Tip Stiffening Cannula available



CM Marked Wire



Detachable Cannula

	SURGIMARC SLIM Qty/Box 10pcs
ĺ	Add <b>X</b> for Double Hook ( VVNSX )
	Add XT for Double Hook and Stiffened Wire ( VVNSXT )
	Add <b>Z</b> for Wire Depth Distal Marking ( VVNSZ )
	Add <b>T</b> for Detachable Cannula ( VVNST )
	Add <b>ZT</b> for Double Hook and Detachable Cannula ( VVNSZT )

CODE						SIZE
<ul><li>VVNS1807</li></ul>	VVNSX1807	VVNSXT1807	VVNSZ1807	VVNST807	VVNSZT1807	18Gx7cm
<ul><li>VVNS1810</li></ul>	VVNSX1810	VVNSXT1810	VVNSZ1810	VVNST1810	VVNSZT1810	18Gx10cm
<ul><li>VVNS2007</li></ul>	VVNSX2007	VVNSXT2007	VVNSZ2007	VVNST2007	VVNSZT2007	20Gx7cm
<ul><li>VVNS2010</li></ul>	VVNSX2010	VVNSXT2010	VVNSZ2010	VVNST2010	VVNSZT2010	20Gx10cm
<ul><li>VVNS2107</li></ul>	VVNSX2107	VVNSXT2107	VVNSZ2107	VVNST2107	VVNSZT2107	21Gx7cm
<ul><li>VVNS2110</li></ul>	VVNSX2110	VVNSXT2110	VVNZS2110	VVNST2110	VVNSZT2110	21Gx10cm

On request Vigeo is available to provide needles with different diameter and length  $% \left( 1\right) =\left( 1\right) \left( 1$ 



One-handed design

## XDRAPE WOMEN'S HEALTH

### Stereotactic Table and Biopsy Station Cover



- Stereotactic Table Drape cover with elastic for quick set-up and removal
- One piece disposable Biopsy Station cover with elastic and adhesive bands
- Fits all the Stereotactic Table and Biopsy Autoguide brands in the market
- Complete fluid containment and easy clean up







XDRAPE		
CODE	DESCRIPTION	Qty/Box
VIG-13-1548EU	Mammotest Table Draping System	20
VIG-10-9601EU	Mammotest Stage Draping System	50
VIG-13-1549EU	Multicare Table Draping System	10
VIG-10-9600EU	Multicare Stage Draping System	20




# SCREEN VISION WOMEN'S HEALTH


### Tissue Transfer and X-Ray System





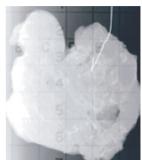

- Safely image and transport of surgical breast specimens without loosing orientation
- Extremely radiolucent and radiopaque grid with pinpointing side indicators
- Precise anatomical assignment of sample margins
- Special stands allowing horizontal and vertical examination
- Once pinned the tumour tissue remains stable and can be easily transported

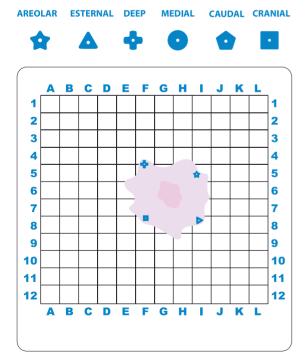


VVX3P Orientation Pins Version

SCREEN VISIO	N	
CODE	THICKNESS	Qty/Box
VVX2	20mm	20
VVX3	30mm	10
VVX3P	30mm + Pins	10



### ISTERIX WOMEN'S HEALTH

### Hysterosalpingography and Hysterosonography Catheter

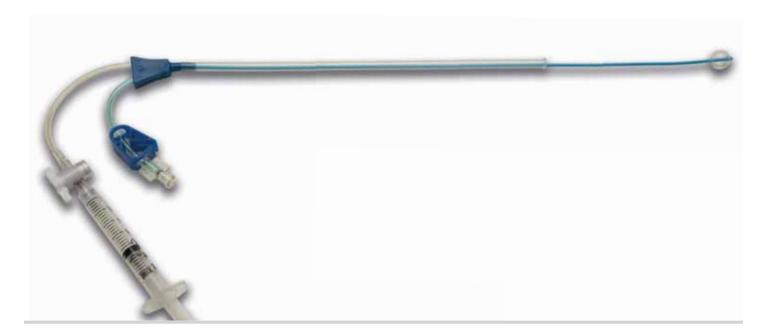

- Flexible radiopaque polyurethane catheter shaft
- Soft, latex-free concentric balloon 3 mm from the catheter tip
- Atraumatic catheter tip
- Front delivery catheter tip version allows guide wire introduction
- Shapeable version with an embedded stylet in cases when stenosis is present or when the vaginal tract and cervical canal are not aligned
- Stable placement sheath for easy introduction of the catheter into the cervical canal
- Syringe included, preset and mounted on the stopcock for balloon inflation port
- Available in two diameters 5F (Balloon 1,5ml) and 7F (Balloon 3,0ml)
- Stylet available for cases when additional stiffness is needed



Shapeable



Front delivery

ISTERIX Qty/Box 10	ocs	
CODE	SIZE	DESCRIPTION
VMI1154	5F	Side delivery
VMI1155	7F	Side delivery
VMI1156	22G	Stylet
VMI1162	5F	Without balloon
VMI1185	5F	Shapeable
VMI1187	7F	Shapeable
VMI1188	5F	Front delivery
VMI1189	7F	Front delivery
VMI2002	5F	Shapeable sheath



# **DUCTAL KIT** WOMEN'S HEALTH

### Galactography Infusion Unit





- Designed to provide ease of insertion with minimal patient discomfort
- Intended for the cannulation of the lactiferous breast duct
- 27G or 30G blunt cannula bent with a 90 degree angle, 300mm PVC tubing and a luer lock hub
- Dilators available ( "D" version )



Mammography

DUCTAL KIT Oty/Box 20pcs			
CODE	SIZE	DESCRIPTION	
VVMK27/C20	27G	Without Dilators	
VVMK30/C20	30G	Without Dilators	
VVMK27/D20	27G	With Dilators	
VVMK30/D20	30G	With Dilators	



### CYTOCAN FINE NEEDLE ASPIRATION

#### Chiba Needle for Aspiration Biopsy and Prenatal Diagnosis



- Puncture cannula developed especially for ultrasonic diagnostics
- Echogenic needle tip assists precise needle tip positioning
- Safe for patient and user, thanks to the adjustable sliding depth stopper
- Graduated puncture cannula to assist in precise depth placement of the needle
- Transparent luer cone for visual inspection of aspirate
- Colour coded hub for easy identification of needle gauge size
- Usable for a variety of applications (cytological biopsy, prenatal puncture, and for percutaneous applications in general)



Handly (code: VVCH)



Syringe Gun (code: SG01)

CYTOCAN Qty/Box 20pcs		
Add <b>U</b> for Chorial kit ( VVBU)		
CODE	SIZE	
<ul><li>VVB1810</li></ul>	18Gx10cm	
● VVB1815	18Gx15cm	
<ul><li>VVB1820</li></ul>	18Gx20cm	
<ul><li>VVB2010</li></ul>	20Gx10cm	
<ul><li>VVB2015</li></ul>	20Gx15cm	
<ul><li>VVB2020</li></ul>	20Gx20cm	
<ul><li>VVB2105</li></ul>	21Gx5cm	
<ul><li>VVB2110</li></ul>	21Gx10cm	
<ul><li>VVB2115</li></ul>	21Gx15cm	
<ul><li>VVB2120</li></ul>	21Gx20cm	
● VVB2210	22Gx10cm	
● VVB2215	22Gx15cm	
● VVB2220	22Gx20cm	
● VVB2310	23Gx10cm	
● VVB2315	23Gx15cm	
<ul><li>VVB2505</li></ul>	25Gx5cm	



## **DPUNCH** SKIN CARE

### Disposable Biopsy Punch



- Sterile single-use, stainless steel instrument for accurate and gentle skin biopsies
- Stainless steel sharp cutting edge
- Ergonomic ribbed handle and anti-slip surface for sure grip control
- Size embossed on each punch for easy identification
- Wide range of sizes



Ergonomic ribbed handle

DPUNCH Qty/Box 20pcs		
CODE	SIZE	
VPT2	2,0mm	
VPT3	3,0mm	
VPT4	4,0mm	
VPT5	5,0mm	
VPT6	6,0mm	



#### Thoracentesis and Paracentesis Set (TAPS)



- Easy,fast, and practical system with single-acting valves
- Waste bags with non-return valve and drainage tap
- Ready to use set with wide choice of needle sizes
- Available in various kit configurations
- 90cm Connection tubes provided with 3-way stopcock
- Disposable, double packaged, sterile



Verres Needle



Needle Catheter

TORAPARAKIT Qty/Box 25pcs				
CODE	DESCRIPTION			
VVV01T1412A	Stopcock 3Way + Verres needle (14Gx12cm) + 1 Drainage bag (2lt)			
VVV01T1412B	Stopcock 3Way + Verres needle (14Gx12cm) + 2 Drainage bag (2lt)			
VVV01T1412C	Stopcock 3Way + Needle Catheter (6Fx12cm) + 1 Drainage bag (2lt)			



# V-DRAIN DRAINAGE

### Single Step Drainage Catheter



- Large oval drainage side holes
- Large lumen and tapered tip
- Locking system available
- Internal diameter of hub and catheter match
- Centimeter depth markings assist in catheter placement
- Echogenic needle tip with Trocar stylet
- Hydrophilic coating for easy penetration
- Easy-to-remove pigtail



Large oval holes



Trocar stylet



Locking system available

CODE	SIZE	
VVVB00620	6Fx19cm	
VVVB00630	6Fx29cm	
VVVB00720	7Fx19cm	
VVVB00730	7Fx29cm	
VVVB00820	8Fx19cm	
VVVB00830	8Fx29cm	
VVVB01020	10Fx19cm	
VVVB01030	10Fx29cm	
VVVB01430	14Fx29cm	
VVVB01630	16Fx29cm	


#### V-DRAIN without locking pigtail Qty/Box 5pcs

CODE	SIZE	
VVVB0620	6Fx19cm	
VVVB0630	6Fx29cm	
VVVB0720	7Fx19cm	
VVVB0730	7Fx29cm	
VVVB0820	8Fx19cm	
VVVB0830	8Fx29cm	
VVVB1020	10Fx19cm	
VVVB1030	10Fx29cm	
VVVB1430	14Fx29cm	
VVVB1630	16Fx29cm	



# VP-JET CLINICAL ONCOLOGY

#### Port Access Needle

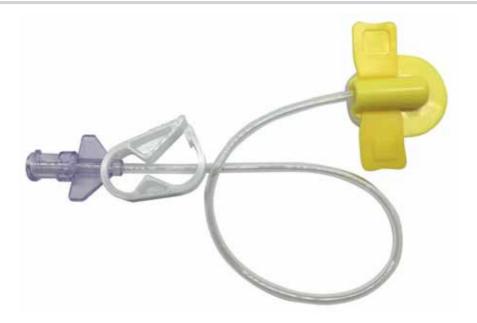

- Single handed infusion port non coring needle access
- Biocompatible materials containing no latex or DEHP
- Retaining wings provide a good grip
- Easy removal of the cannula
- Luer lock Y-site allows additional medication
- Cushioned platform for patient comfort





Y-site allows additional medication



Cushioned platform

VP-JET Qty/Box 50µ	ocs
Add <b>Y</b> for double	le injection site
Add <b>YS</b> for dou	ble injection site needle free
CODE	SIZE
<ul><li>VP2020</li></ul>	20Gx20cm
<ul><li>VP2025</li></ul>	20Gx25cm
● VP2220	22Gx20mm
● VP2225	22Gx25cm



# ALCOJECT CLINICAL ONCOLOGY

Percutaneous Ethanol Injection Needle (PEIT)



- Larger volume of treated tissue
- Lower ethanol flow into the peritoneal cavity
- No defelction of the needle from the path
- Suitable for small and deep nodules
- Centimeter depth markings assist in needle placement
- Indicator shows flow direction
- Colour coded for quick gauge size identification
- Luer extension tube for a flexible connection to the syringe



On request Vigeo is	s available to	provide	needles	with	different	diameter	and I	ength

• VVHA32115 21Gx15cm Open Tip

SIZE

ALCOJECT Qty/Box 10pcs

CODE



DESCRIPTION

# **BLOCK** CLINICAL ONCOLOGY

### Bone and Bone Marrow Biopsy Device



- Designed for quick and easy cortical penetration
- Ergonomic T- handle for increased control
- Stylet with pyramidal tip
- Cannula with fish mouth tip
- Centimeter depth markings assist in needle placement
- Extraction cannula with centimeter markings and blunt tip
- Improved patient comfort
- B-Locking technology to avoid deflection during the sample release



Fish mouth tip



B-Locking technology

BLOCK Qty/Box 10pcs				
CODE	SIZE			
● VVJT0810	8Gx10cm			
● VVJT0815	8Gx15cm			
○ VVJT1110	11Gx10cm			
OVVJT1115	11Gx15cm			
<ul><li>VVJT1307</li></ul>	13Gx7cm			
<ul><li>VVJT1310</li></ul>	13Gx10cm			



# **BCUT** CLINICAL ONCOLOGY

### Bone and Bone Marrow Biopsy Device



- Designed for quick and easy cortical penetration
- Ergonomic T- handle for increased control
- Stylet with pyramidal tip
- Cannula with fish mouth tip
- Centimeter depth markings assist in needle placement
- Extraction cannula with centimeter markings and blunt tip
- Improved patient comfort



Fish mouth tip

BCUT Qty/Box 10pcs			
CODE	SIZE		
<ul><li>VVJ0810</li></ul>	8Gx10cm		
<ul><li>VVJ0815</li></ul>	8Gx15cm		
○ VVJ1110	11Gx10cm		
○ VVJ1115	11Gx15cm		
<ul><li>VVJ1307</li></ul>	13Gx7cm		
<ul><li>VVJ1310</li></ul>	13Gx10cm		



# MYELO CLINICAL ONCOLOGY

### Bone and Bone Marrow Biopsy Device



- Cannula and stylet designed to allow an easier penetration of the bone cortex
- Adjustable depth stop for sternal or paediatric application, which can be removed for iliac crest use
- Luer slip and luer-lock syringes fit securely



MYELO Qty/Box 10pc	cs
Add 1 for Ergono	mic handle ( VVI1 )
Add <b>E</b> for Explan	t version ( VVIE )
CODE	SIZE
VVI1410/50	14Gx10/50mm
VVI1430/70	14Gx30/70mm
VVI1510/50	15Gx10/50mm
VVI1530/70	15Gx30/70mm
VVI1610/50	16Gx10/50mm
VVI1630/70	16Gx30/70mm
VVI1810/50	18Gx10/50mm



# ERGO CLINICAL ONCOLOGY

### Bone and Bone Marrow Biopsy Device



- Ergonomic T-handle design to minimize hand and arm fatigue
- Cannula and stylet designed to allow an easier penetration of the bone cortex
- Adjustable depth stop for sternal or paediatric application, which can be removed for iliac crest use
- Luer slip and luer-lock syringes fit securely



Colour code

ERGO Qty/Box 10pcs	
Add <b>E</b> for Explant	version ( VVILE )
CODE	SIZE
VVIL1410/50	14Gx10/50mm
VVIL1430/70	14Gx30/70mm
VVIL1510/50	15Gx10/50mm
VVIL1530/70	15Gx30/70mm
VVIL1610/50	16Gx10/50mm
VVIL1810/50	18Gx10/50mm



## T'AM KIT CLINICAL ONCOLOGY

### CT-Guided Percutaneous Biopsy of Spinal Lesions





- Chiba fine needles for local anesthesia included
- Seldinger placement of introducer
- Suitable for accurate biopsy of deep bone lesions
- An inner trocar guides the outer cannula to the lesion margin offering control and precision whilst traversing dense bone
- The cannula remains in situ permitting multiple biopsies
- High grade steel and a Trephine tip, facilitate superior sampling of critical lesions



CT

T'AM KIT Qty/Box 5 <i>pcs</i>				
CODE	SIZE			
TCAK2510	13Gx100mm			
TCAK2512	13Gx125mm			
TCAK2515	13Gx150mm			
TCAK3010	11Gx100mm			
TCAK3125	11Gx125mm			
TCAK3015	11Gx150mm			



## BIOMARC IGRT RADIOTHERAPY

### Exclusive Distributor C 6 0086








- Used to target localization for lesions or organs which move in respect to external anatomy
- Pure gold cylinders version specially knurled to inhibit migration
- Pyrolytic carbon coated zirconium oxide bar-bell version to reduce the artifacts
- Provide visualization under fluoroscopy, X-ray, CT, MRI, ultrasound and other imaging modalities
- Preloaded in a pre-waxed tipped needle or in a wax- free special version
- Each individually packaged for wide flexibility during the treatment



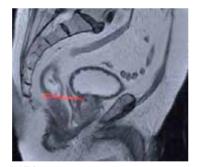

CT MRI

BIOMARC IGRT			
CODE	SIZE	DESCRIPTION	Qty/Box
040116	17Gx20cm	1,2x3mm Gold Marker	3
040118	17Gx20cm	1x4mm Gold Marker-Ergonomic handle	10
040125	17Gx30cm	1x3mm Gold Marker	3
040117	18Gx20cm	1x5mm Ceramic Marker	1
040135	18,5Gx13cm	1x3mm Ceramic Marker	3
040113	18Gx20cm	1x3mm Ceramic Marker	3
040146	18Gx20cm	0,9x5mm Ceramic Marker- No Artifact	3



# ARCTIC - COLDSPOT INFECTION CONTROL

### Sterile Cold Therapy





- Available Pre-Frozen or Instant Cold pack
- Reduces bleeding, bruising, pain, and inflammation
- Suitable after minimally invasive Breast and Soft Tissue Biopsy procedures



After breast biopsy

ARCTIC - COLDSPOT			
CODE	DESCRIPTION	Qty/Box	
500	COLD SPOT Pre-frozen Gel	100	
V500	ARCTIC Instant Cold Pack	40	



### **ECOVER** INFECTION CONTROL

#### **Transducer Covers**











Disposable covers for US procedures

- General purpose sterile set includes one PU probe cover, one Sterile gel packet and two elastic bands.
- Premium version provides quick and easy loading of an ultrasound probe into a long, sterile probe cover
- Safe adhesive version eliminates the need of Gel under the Cover



General Purpose set



Premium Version



Safe Adhesive version

ECOVER Qty/Box 20pcs
----------------------

Add **21 or 22** for General Purpose set (VIG-11-21 / VIG-11-22)

Add 12 for Premium Version (VIG-11-12)

Add 11 for Safe Adhesive Version (VIG-11-11)

Add II for Sale Adriesive version ( vig-11-11 )		
CODE	SIZE	DESCRIPTION
VIG-11-2102EU	10x61cm	PU Cover + Sterile Gel + Elastic Band
VIG-11-2105EU	15x122cm	PU Cover + Sterile Gel + Elastic Band
VIG-11-2107EU	15x244cm	PU Cover + Sterile Gel + Elastic Band
VIG-11-2251EU	14x61cm	PU Cover + Sterile Gel + Elastic Band
VIG-11-2253EU	14x91cm	PU Cover + Sterile Gel + Elastic Band
VIG-11-2252EU	14x147cm	PU Cover + Sterile Gel + Elastic Band
VIG-11-1251EU	14x61cm	PU Cover Kit with ergonomic handle
VIG-11-1253EU	14x122cm	PU Cover Kit with ergonomic handle
VIG-11-1139EU	18x127cm	PE Cover with Adhesive
VIG-11-1140EU	18x249cm	PE Cover with Adhesive



## **DECOMAT** INFECTION CONTROL

### Multilayer Decontamination Mat



- Used in all environments that require the total absence of contaminating agents and the preservation of high hygiene levels
- Barrier against the growth and spreading of microorganism
- Antibacterial and antifungal glue captures and retains dust and bacteria
- 2mm thickness Mat includes 30 adhesive PE sheets
- Non-adhesive strip for easy removal of used layers
- Dispenser folder for confortable handling and vertical storage




<b>DECOMAT</b> Qty/B	ox 5 <i>pcs</i>
CODE	SIZE
VK450	45x120cm
VK600	60x120cm
VK900	90x120cm





₩ **-**

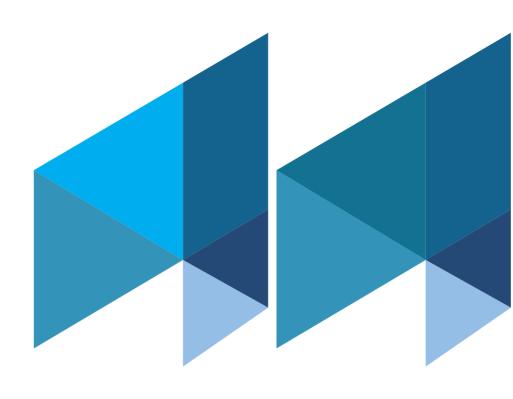

UNI CEI EN ISO 13485:2012

May 2018

www.vigeosrl.it





# vigeo

VIGEO srl
San Biagio - Via Dell'Alpino 2
46031 Bagnolo San Vito MN - ITALIA
Tel. +39 0376 449097 - Fax +39 0376 282245
www.vigeosrl.it
customerservice@vigeosrl.it